
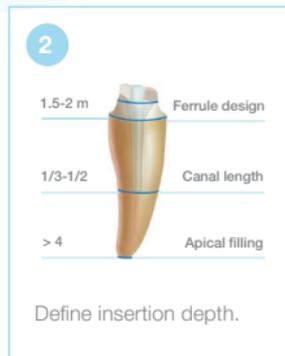


## D.T. Light-Posts – specific Double Taper Shape

Step by Step  
Card

# Placing root canal posts with the D.T. Light-Post System\*

## A. Preparation of Post Space and D.T. Light-Post

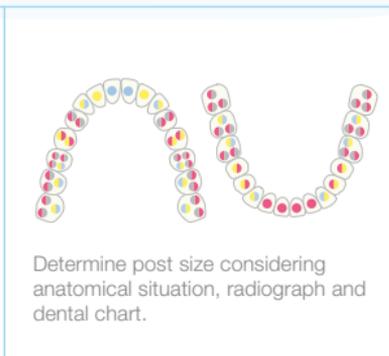

## B. Adhesive Bonding

8

**optional**



Etch, rinse and dry both canal and coronal dentin.\*

9



Recommended product:  
Prime&Bond active®  
and Self Cure Activator  
from Dentsply  
Sirona

Mix adhesive and Self Cure Activator, apply mixture and leave undisturbed for min. 20 sec. Remove excess. Gently air dry.\*

10



Lightly coat with adhesive/SCA mixture and gently air dry.

11



Recommended product:  
core-X® flow  
from  
Dentsply Sirona

Apply composite cement into the root canal.\*

12



Place D.T. Light-Post, stabilize and light cure.\*

13



Create core build up with composite cement.\* Cover post completely.

\*Please read and follow the instructions for use of D.T. Light-Posts and adhesive luting agents used



Endo  
Easy  
Efficient

For further information please contact your preferred distributor or visit [www.vdw-dental.com](http://www.vdw-dental.com).

 R.T.D. (RECHERCHES TECHNIQUES DENTAIRE)  
3 Rue Louis Néel, Technoparc, Espace Gavanière  
38120 Saint-Égrève  
France  
Tel +33 4 76 56 56 66  
Fax +33 4 76 56 11 77  
[info@rtd.fr](mailto:info@rtd.fr)

Distributed by:

**VDW GmbH**  
Bayerwaldstr. 15  
81737 Munich  
Germany  
Tel +49 (0) 89 62734-0  
Fax +49 (0) 89 62734-304  
[info@vdw-dental.com](mailto:info@vdw-dental.com)

  
0459