

BeeFill™ 2in1

Time-saving, clean and long-lasting 3-dimensional obturation of root canals

Best obturation results for every tooth

- Best filling results with warm, vertical compaction (Toronto Study, 2003 & 2004)
- Vertical compaction is a faster method to achieve a long-lasting, tightly sealed, 3-dimensional root canal filling
- Reliable obturation of lateral canals
- Minimal risk of root fractures
- Guttapercha filling to desired level for immediate placement of a fiber post
- Adhesive restoration in one appointment



→ Filling was never so easy

- Downpack and Backfill function in one compact device
- Pulp sensitivity testing
- Simple navigation at the push of a button
- All parameters individually adjustable
- Clear display and controls for Downpack or Backfill
- Two preset buttons allow preferred settings to be saved




→ Designed for professional obturation


- Ergonomically designed handpiece with a 360° operation angle
- Feather-touch activation by flexible silicone ring
- Excellent sight of the root canal thanks to slim handpieces and extra-long obturation cannulas
- Micromotor controls precise release of guttapercha
- Visual and audible signals for operation and safety on the handpieces
- Blue LED lights indicate active handpiece




→ Comfortable filling

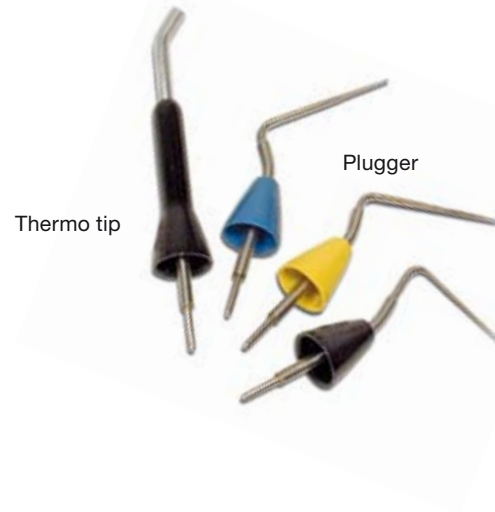

- Quick change of cartridges through drop-in/drop-out design
- Extra-long, flexible cannulas, suitable for curved root canals
- One cartridge fills 4-6 root canals
- Standard and soft pluggers
- ISO coloured pluggers with length markings for downpacking or searing off Guttapercha
- Thermal response tip for pulp sensitivity testing

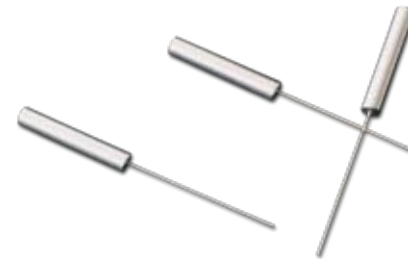

Guttapercha cartridges

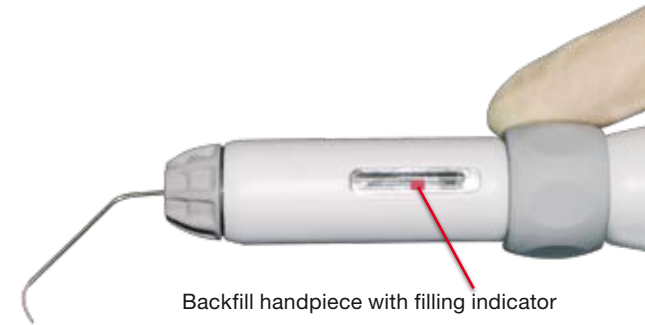

Range of products



Ref. 0541 BeeFill™2in1 obturation device

Ref. 0545 BeeFill™2in1 obturation kit

- 1 obturation device
- 1 plugger small ISO 40/.03
- 1 root canal sealer 2Seal
- 10 guttapercha-cartridges 20G
- 10 guttapercha-cartridges 23G
- 2 Machtou-Plugger, sizes 1-2 and Gr. 3-4

BeeFill™ Plugger standard

Ref. 0547 000 040 ISO 40 Taper .03 (1 piece)

BeeFill™ Plugger soft

Ref. 0546 000 050 ISO 50 Taper .05 (1 piece)

Ref. 0546 000 060 ISO 60 Taper .06 (1 piece)

BeeFill™ Thermo tip

Ref. 0548 000 000 Pulp sensitivity tester (1 piece)

BeeFill™ Guttapercha Cartridges

Ref. 0532 053 080 20G / 0.8mm (10 pieces)

Ref. 0532 053 060 23G / 0.6mm (10 pieces)

Ref. 0531 BeeFill™ Backfill device

Ref. 0535 BeeFill™ Backfill Kit

- 1 Backfill device
- 10 cartridges 20G
- 10 cartridges 23G
- 2 Machtou-Plugger, sizes 1-2 and 3-4

Ref. 0544 BeeFill™Pack Downpack device

Ref. Nr. 0544 000 600

BeeFill™Pack Downpack Kit

- 1 BeeFill™ Pack
- 1 Plugger small ISO 40 Taper .03
- 1 root canal sealer 2Seal

